

FREELIFT Classic seat

Single tube curved stairlift system



Main Benefits

If you are thinking about installing a stairlift, and your stairs turn corners or are spiral, Handicare's Classic seat and single tube stairlift system offers you a safe and reliable solution.

Handicare's stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.

Handicare's single tube stairlift system suits most staircases with bends or corners and when installed with the Classic seat it offers you all the basic functions that you would expect plus the option to upgrade to a powered swivel and/or have the fold up footplate linked to the seat so that folding the seat will automatically fold the footplate up too.

With a choice of three leather look upholstery colours and three standard track colours the Classic seat single tube stairlift lift system will compliment


schemes.

All Handicare stairlifts have a battery backup which means that your stairlift will continue to work in the event of power cut so you do not need to worry about being stranded up or downstairs.

The compact unit folds up leaving the staircase free for other users and the bespoke track is designed to fit your individual staircase and maximise clear space on the stairway. It can be supplied with a hinged track if you need to avoid obstructing doorways close to the bottom of your stairs.

Fitted to the stairs, not the wall, it is unlikely you will need to redecorate following the installation of your Handicare stairlift.

Your Handicare stairlift will also come with two remote controls which means that, once folded, you can easily send it to the upstairs charge point freeing up space in your hallway.

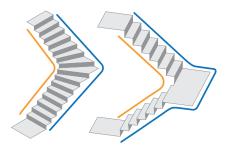

Key: - Rembrandt -

- Van Gogh





Which track?

Every staircase that is curved or turns a corner will have a shorter and a longer side. You may hear stairlift tracks being referred to as inside and outside or internal and external depending on which side to the staircase they sit on.

Handicare's single tube system has two names: the 'Rembrandt' travels along the shorter side of your staircase (the internal) and the 'Van Gogh' travels the longer side.

Your stairlift supplier will advise on which is the best solution for your staircase.


Safety as standard


Your Handicare stairlift will come with safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up the stairs and will come with a retractable seatbelt as standard.


Easy to use

The operating joystick is easy to use, simply hold it in the direction you wish to travel, up or down, and your Handicare stairlift will move smoothly and gently on its way.



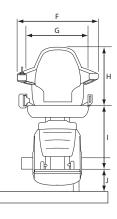


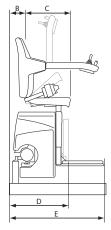

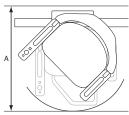

Classic	Van Gogh	Rembrandt
A Swivel radius*	645-715mm	665-735mm
B Minimum wall to seat back	100mm	120mm
C Seat depth	370mm	370mm
D Minimum folded width to footplate	404mm	424mm
E Minimum open width to edge of footplate	624mm	644mm
F Armrest width – external	555mm	555mm
G Armrest width – internal	439mm	439mm
H Seat back height	430mm	430mm
I Footplate to seat height	452mm	452mm
J Minimum footplate height	133mm	133mm
Bottom parking space requir ed	650mm	650mm
Minimum track intrusion into staircase	170mm	190mm
Maximum weight limit	125kg	125kg
Speed	0.11m/sec (approx)	
Power supply	24v DC battery	
Charger	100-240v (50-60Hz)	
Maximum staircase angle	60°	

^{*}Dimensions vary depending on hand of seat

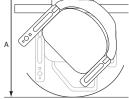







RED















Norfolk Stairlifts, Unit 10, Ironside Way, Hingham, Norwich, Norfolk NR9 4LF.



Tel: 01953 850004

Norfolk Stairlifts is an approved distributor for this manufacturer

Please call 01953 850 004 for more information or a no-obligation FREE survey.

- Established in 1998 as the only independent specialist in Norfolk.
- Trained to install most makes /models, so you can choose what suits you.
- Qualified to repair all makes of stairlifts should you need us.
- Passionate about serving our 700 satisfied customers in Norfolk & Suffolk.

*Please check if this particular model is on display before arriving at the showroom.

